



A refined design. An economical choice.

InnoSpire Essence compressor nebulizer system

When cost is a primary consideration, choose InnoSpire Essence. Essence, part of the InnoSpire family of compressor nebulizer systems from Philips Respironics, provides fast and efficient aerosol drug delivery in a compressor nebulizer system that's priced right.

Combined with proven SideStream technology, Essence is compact, lightweight, and economical. Its stylish design incorporates a handset docking station, filter and power cord and blends gracefully into the home environment. It is covered by a 5-year warranty and is compliant with the international safety and essential performance standard IEC 60601-1 third edition and its amendments.

PHILIPS

RESPIRONICS

sense **and** simplicity


InnoSpire Essence



For custom logo compressors contact your Philips Respironics sales representative.

SideStream disposable and reusable nebulizers

Essence works efficiently with SideStream nebulizers that feature an active venturi system and provide air flow in addition to that supplied by the compressor, resulting in faster drug delivery and shorter treatment times. They are easy to use with fewer parts to assemble or clean than other conventional nebulizers. SideStream's unique design provides consistent delivery of medication time after time.

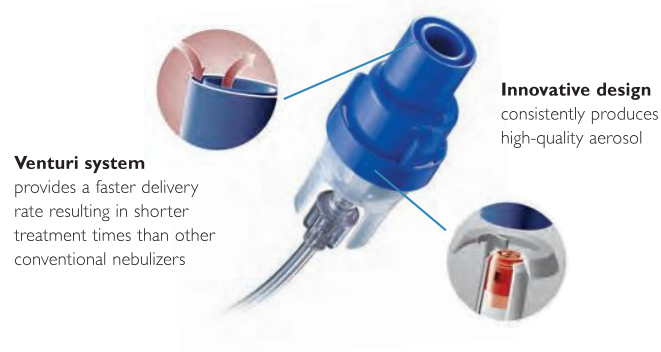

Product specifications

Weight	3.3 lbs
Size	4.2" x 6.5" x 6.5"
Electrical rating/power	120V/60Hz, 1.7 Amp
Nebulization capacity	8 ml
MMAD	(80%) <5 microns
Treatment time	5-7 minutes
Maximum pressure	44 psi
Maximum flow rate	9.4 LPM
Operating flow rate	7 LPM @ 10 psi
Mode of operation	Intermittent (30 min. on /30 min. off)
Warranty	5 years

Ordering information

Product description	Reorder number	NDC number
InnoSpire Essence with SideStream disposable nebulizer	1099966	08373-9966-00
InnoSpire Essence with SideStream disposable and reusable nebulizers	1100312	08373-0312-00
InnoSpire Essence with SideStream disposable and reusable nebulizers, Tucker the Turtle mask and carrying case	1100313	08373-0313-00
Replacement filter (4-pack)	1102088	08373-2088-00
Standard carrying case	1101452	08373-0452-00
Premium carrying case	1026568	08373-0568-00

Economical

Product description	HCPCS
InnoSpire Essence compressor nebulizer system	E0570
SideStream disposable nebulizer	A7003
SideStream reusable nebulizer	A7005

Please visit www.philips.us/essence for more information

Respironics, InnoSpire, Essence, SideStream and Tucker the Turtle are trademarks of Koninklijke Philips Electronics N.V. and its affiliates. All rights reserved.

© 2012 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication. CAUTION: US federal law restricts this device to sale by or on the order of a physician.

DD MCI 4104968 PN 1101687

Philips Healthcare is part of Royal Philips Electronics

Philips Respironics International
Headquarters
+33 1 47 28 30 82
Philips Respironics Asia Pacific
+65 6882 5282

Respironics New Jersey, Inc.
5 Wood Hollow Road
Parsippany, NJ 07054
+1 800 345 6443
+1 724 387 4000
www.philips.com/respironics

